
User Manual

Naipo Cordless Handheld Percussion Massager

Shenzhen Pinlong E-Commerce Co., Ltd.
 Web: www.naipocare.com
 E-mail: support@naipocare.com
 Add: Huanan City Electronic Trading Building P05
 Room 116, Longgang District, Shenzhen, China.

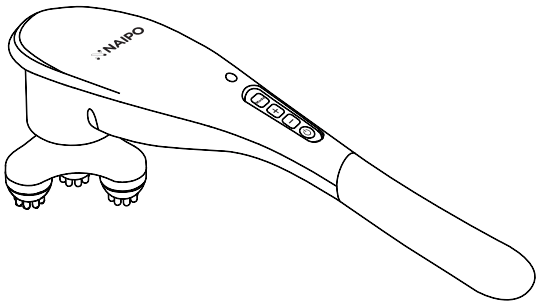Model: MGPC-5610

Contents

English ----- 01-04
 Deutsch ----- 05-08
 Français ----- 09-12
 Español ----- 13-16
 Italiano ----- 17-20

Introduction

Thank you for purchasing the Naipo Cordless Handheld Percussion Massager. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.

Package Included

- 1 x Handheld Massager
- 1 x User Manual
- 3 x Interchangeable Attachments
- 1 x Charger
- 1 x Stand

Product Specs

Model	MGPC-5610
Name	Naipo Cordless Handheld Percussion Massager
Power input	DC 9V
Power consumption	9 W
Default time	20 mins
Li-ion	1800 mAh
Charging Time	2.5 h
Usable Time	1.5 h

Warnings

- Please pull out the plug from the socket when cleaning, maintenance or moving.
- Please turn off the device if there is something unusual occurs.
- Please do not use wet hand to pull out or insert the power cord.
- Please do not use it in humid places such as the bathroom.


- 01 -

- Please do not destroy, bend or handle the power line, or it may cause fire.
- Please clear the dust of adaptor regularly.
- Please avoid expose to sunlight for long time.
- Please do not modify the product casually.
- Please do not use it while sleeping.
- Please do not wash the product by water.
- The product should be supervised by people if used for long time.
- Children are prohibited from using this product, please put it in a safe place which is out of their reach.

Q&A


Problem	Remedy
The device doesn't work.	Please check if the power supply is normal. Please reconnect the plug to the mains supply.
During operation, the controller is working while the massager is not.	The massager is overheated because of long time use. Please disconnect device from mains power, wait at least 50 mins before switching on. Press the power button to restart the device.
The device stops working suddenly.	Please check if the power supply is normal. Please reconnect the plug to the mains supply.

Product Diagram

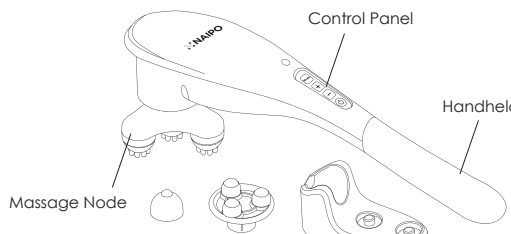- 02 -

Controller Diagram

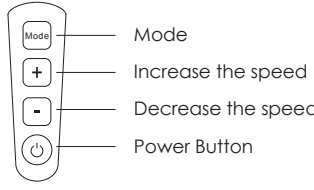Product Operation Guide

Operation preparation

- Select the proper massage node.
- Connect the adapter correctly and plug it into the power supply.

Instructions for use

Power Button

Press the power button to activate the vibration function, press it again to turn off the device.

[+]
Press the [+] button to increase the vibration speed, 5 levels of speed are available.

[-]
Press [-] button to decrease the vibration speed.

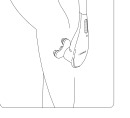
Mode Button

Press mode button to activate auto massage program, 5 auto programs are available (change from compact to slow).

- 03 -

Massage Areas

Warranty & Customer Support

For questions, assistance, support, or warranty claims, please contact us at the address below that corresponds with your region. Please include your Amazon order number and product model number.

Website: www.naipocare.com

Email: support@naipocare.com

*Please note: Naipo can only provide after sale service for products purchased directly from Naipo. If you have purchased from a different seller, please contact them directly for service or warranty issues.

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