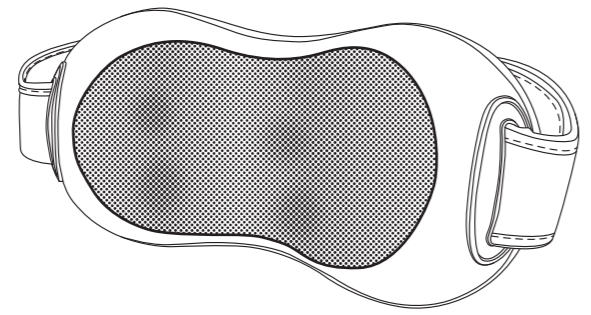
Shenzhen Pisheng E-Commerce Co., Ltd.
Web: www.naipo.com
E-mail: support@naipo.com
Add: Huanan City Electronic Trading Building P05
Room 116, Longgang District, Shenzhen, China.

CE **RoHS** 

NAIPO
SINCE 2005

User Manual

Shiatsu Pillow with Kneading and Heat

Model: MGP-999

Contents

English	01-04
Deutsch	05-08
Français	09-12
Español	13-16
Italiano	17-20

Introduction

Thank you for purchasing the Naipo Shiatsu Pillow with Kneading and Heat. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.

Package Included

- 1 x Shiatsu Pillow
- 1 x User Manual
- 1 x Car Charger
- 1 x Charger

Product Specs

Model	MGP-999
Name	Naipo Shiatsu Pillow with Kneading and Heat
Power Input	AC 100V-240V / DC 12V
Power consumption	24 W
Rated frequency	50 / 60 Hz
Default time	15 mins
Security Level	II

Important Notes

If you are suffering from one of the health complications or conditions below, please consult your doctor before use:

- If you have cancer.
- If you have any blood-related problems.
- If you have heart disease.
- If you have any spinal problems.
- If you are in poor physical health.
- If you feel sick or nauseous.
- If you are pregnant or menstruating.
- If you have pacemaker or other medical implants.


Children are prohibited from using this product; please put it in a safe place which out of their reach.


People with physical, sensory or mental disabilities should not use the product.

QLA

Problem	Remedy
The device doesn't work.	Restart the massage after it cools down.
The intensity of massage is too strong or too weak while using.	Adjust the angle until the intensity is moderate. Adjust the intensity that you lean on the massage.

Product Diagram

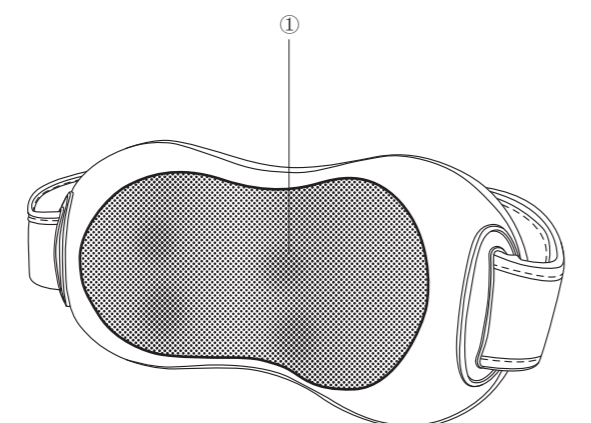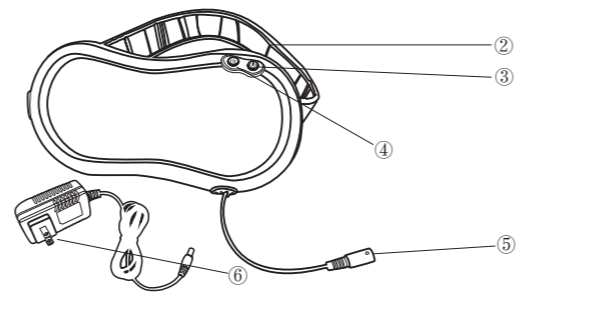
① Massage Node
② Headage
③ Power Button
④ Head Button
⑤ Power Connector
⑥ Power Plug

Product Operation Guide

- To use the product in your home, connect the charger and the product, then plug it into the main power outlet.
- To use the product in your car, connect the car charger and the product, then plug it into the cigarette lighter socket.
- Press the power button once, the product will activate auto massage (kneading and heat). And the kneading function will change mode for every 1 minute.
- Press the mode button to change another massage mode.
- Press the heat button once and the massage speed will become quicker. Press it again to revert to the original state.
- Long press the heat button for 2 seconds to turn off the heat function, and long press it again to turn on the heat function.
- Long press the power button for 3 seconds to turn off the device.
- The product will turn off automatically after 15 minutes. We recommend you to turn it off at least least 30 minutes after each use.

Massage Areas

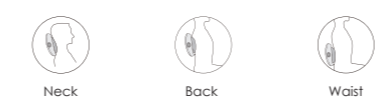Neck Back Waist

Warranty & Customer Support

For questions, assistance, support or warranty claims, please contact us at the address below that corresponds with your region. Please include your Amazon order number and product model number.

Address: Huanan City Electronic Trading Building P05 Room 116, Longgang District, Shenzhen, China.
Website: www.naipo.com
E-mail: support@naipo.com

*Please note: Naipo can only provide after sale service for products purchased directly from Naipo. If you have purchased from a different seller, please contact them directly for service or warranty issues.

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