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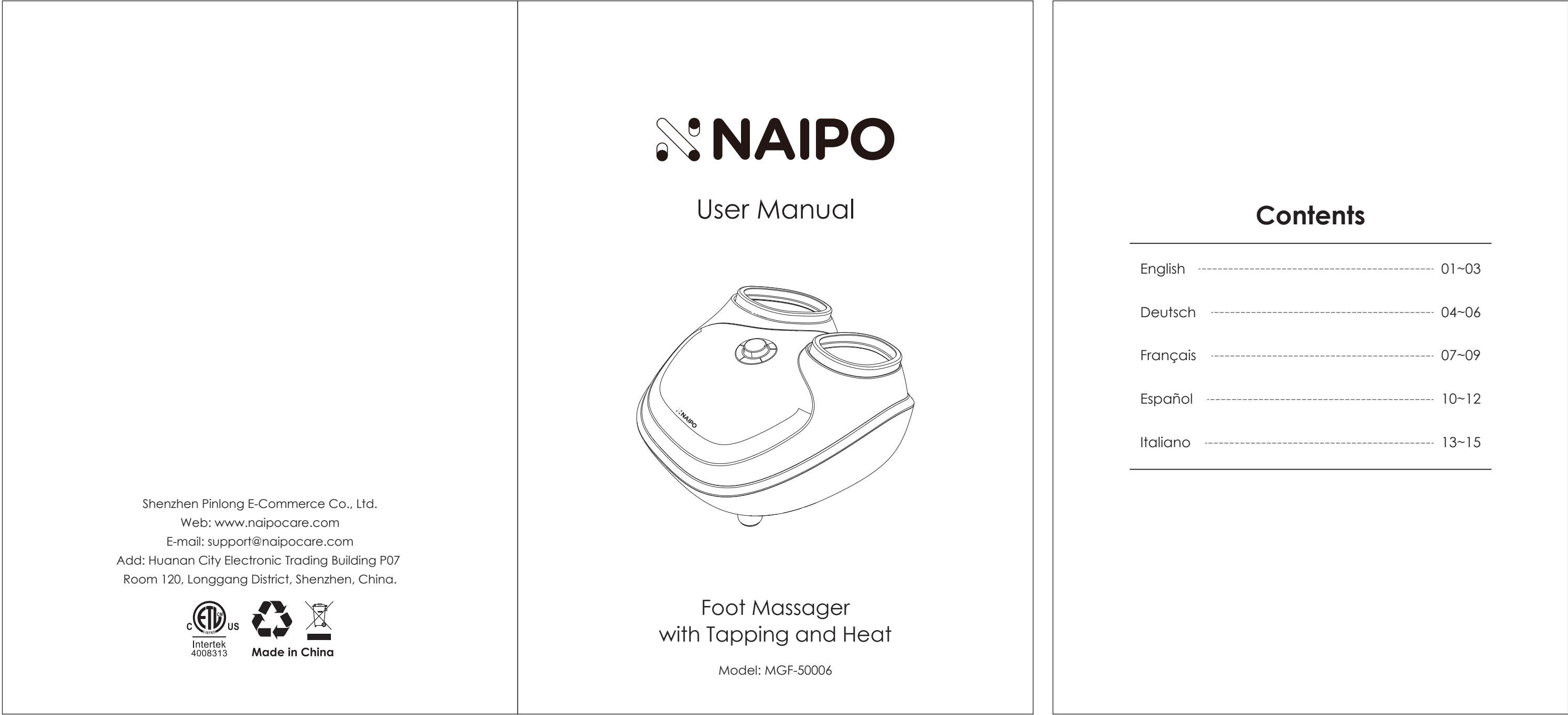
Introduction

Thank you for purchasing the Naipo Foot Massager with Tapping and Heat. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.

Included

1x Foot Massager
1x User Manual

Product Specifications


Product Name	Naipo Foot Massager with Tapping and Heat
Model	MGF- 50006
Power Source	120 V~
Power Frequency	60 Hz
Default Timer	15 Minutes
Dimensions (Length/Width/Height)	410 x 390 x 261 (mm) 16.1 x 15.4 x 10.3 (in)
Weight	7.65 KG / 16.87 Pounds

Important Safeguards

- Do not use the product in a wet or dusty environment as this may result in malfunction or electric shock.
- Do not put the product in bed.
- Do not place the product near liquids.
- Do not sit on the product.
- Do not expose the product to direct sunlight.
- Do not use the product if you are experiencing the following situations: heart attacks, pregnancy, or any other health concerns.
- Use of this product for more than 15 minutes per time is not recommended.
- The product is prohibited from being used by children without supervision from adults.
- Do not use this product if you are going to sleep soon.
- Do not use the product if the foot cloth is torn out.
- Do not use the product if water permeates the product.
- Stop using the product if abnormal or severe pain occur.
- Do not use the product if the product works abnormally.


- 01 -

Product Diagram

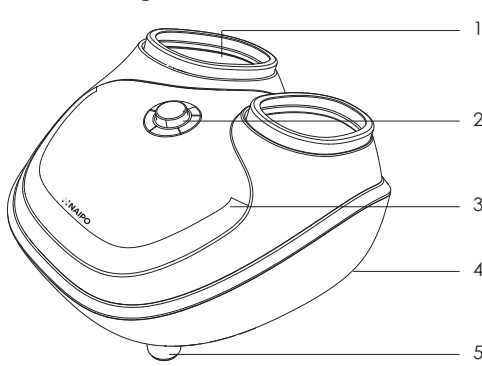1. Washable foot cloth
2. Control panel
3. Main body
4. Cable storage set
5. Footstand

Control Panel

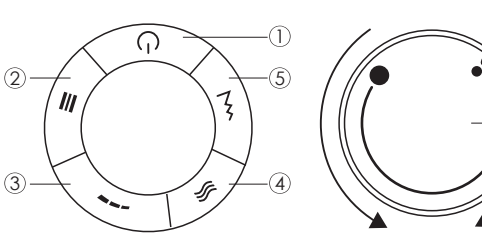- ① Power
- ② Mode
- ③ Air compression intensity
- ④ Heat
- ⑤ Tapping
- ⑥ Tapping speed control

Product Operation Guide

Note: Before being put into use, make sure that the product is in good condition without any damages and that it's positioned flat. Connect the product to the power supply and there will be an audible signal. After the product has been turned on, there will be audible signals for the press of the buttons each time.

1) Power

Press the button once to revert the product to standby mode (red light) with white round LED light on and twice to turn off it (The white light is always on until the product is powered off).

2) Mode - 3 default modes

- a) Mode 1 (green light): air compression (weak) + heating + tapping (intermittent) + rolling (intermittent)
- b) Mode 2 (blue light): air compression + heating + tapping + rolling
- c) Mode 3 (red light): air compression + heating + rolling+ tapping (intermittent)

3) Air Compression Intensity

When the default mode is activated, the air compression intensity is low by default (green light). Press the button once to activate medium-pressure air compression (blue light) and twice to activate high-pressure air compression (red light).

- 02 -

4) Heat

- a) When tapping function is activated, press the button once to activate heating function (red light) and twice to turn off the function.
- b) When default mode is activated, the heating function is on by default. Press the button once to turn off heating function.

5) Tapping - 3 alternative modes (can be activated separately under standby mode)

Under standby mode, press the button once to activate consistent tapping massage (green light), twice to activate intermittent tapping massage (blue light), and three times to activate multi-frequency tapping massage (red light).

6) Tapping Speed Control

Rotate clockwise to increase the speed and counterclockwise to decrease the speed.

7) Timer

When reverting to standby mode, the product will automatically shut off after 15 minutes of use. Disconnect from the power supply when the product is not in use.

Product Maintenance

- Disconnect the product from the power supply before cleaning.
- Clean the product only with a soft and slightly damp sponge.
- Avoid cleaning the product with any corrosive detergent, such as Benzene or diluents.

Q&A

Problem	Remedy
The product doesn't work.	1. Please connect the product to the power supply. 2. Please press the power button. 3. Remove the protective film of the control panel.
The product stops working while in use.	1. Please make sure the product is connected to the power supply correctly. 2. Press the power button again after the product has cooled down. 3. The product has been in consistent use for more than 15 minutes.

Warranty and Customer Support

For questions, assistance, support, or warranty claims, please contact us at the address below that corresponds with your region. Please include your Amazon order number and product model number.

Address: Huanan City Electronic Trading Building P07 Room 120, Longgang District, Shenzhen, China.
Website: www.naipocare.com
E-mail: support@naipocare.com

*Please note: Naipo can provide after sale service for products purchased directly from Naipo. If you have purchased from a different seller, please contact them directly for service or warranty issues.

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