


User Manual

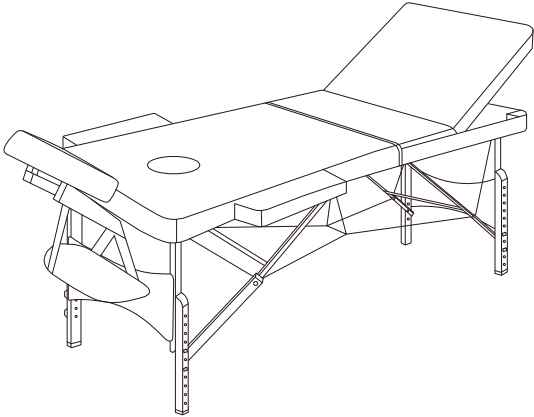

Three-section Massage Table with Wooden Legs

Model: MGBC-1302

Introduction

Thank you for purchasing the Naipo Three-section Massage Table with Wooden Legs. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.


Included

- 1 x Massage table
- 1 x Face hole plug
- 1 x Face cradle cushion
- 1 x Face cradle
- 1 x Armrest
- 2 x Side armrest
- 1 x Carrying bag
- 1 x User manual

Product Specifications

Model	MGBC-1302
Name	Naipo Three-section Massage Table with Wooden Legs
Dimensions	185*60*(61~88)cm (adjustable height)
Working weight capacity	Less than 225kg

Using Guide

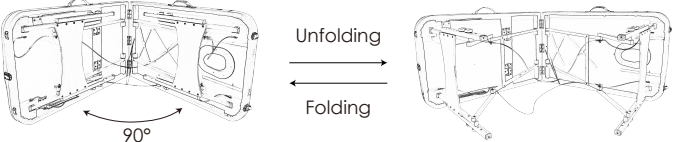


- EN -

User Tips:

- a. When unfolding, lay the folded product on its side, unfold it until the angle is 90°, and then unfold the wooden legs with two hands.
- b. When folding, lay the unfolded product on its side, fold the wooden legs, and then fold the product slowly.
- c. Please make sure that the wire ropes do not wind around the wooden legs when unfolding or folding the product.

Installation of Face Cradle Cushion

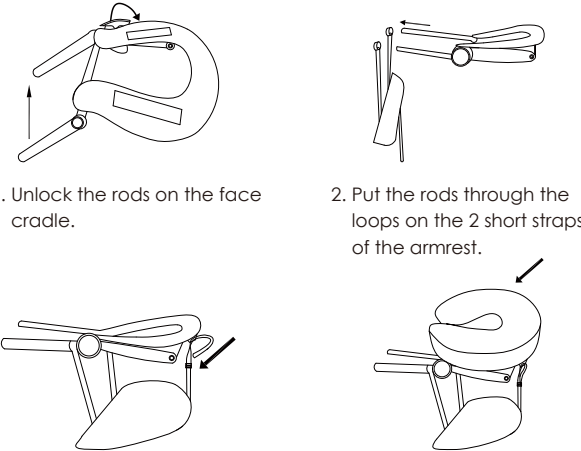

1. Unlock the rods on the face cradle.





2. Put the rods through the loops on the 2 short straps of the armrest.

3. Put the long strap of the armrest through the hoop of the face cradle and stabilize with the back buckle.


4. Put the face cradle cushion on the face cradle.

- EN -

Shenzhen Pinlong E-Commerce Co., Ltd.
Web: www.naipocare.com
E-mail: support@naipocare.com
Add: Huanan City Electronic Trading Building P07
Room 120, Longgang District, Shenzhen, China.




Height Adjustment

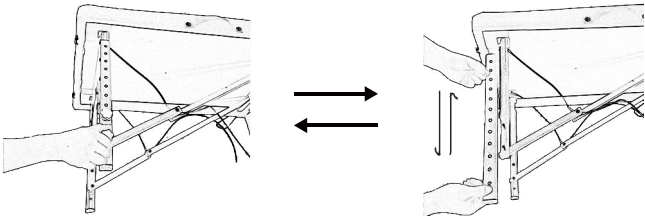

Rotate the 2 knobs counterclockwise on a wooden leg until they fall down. Take down the wooden leg and put it on the desired height with the screws through the corresponding holes. Then put the knobs back and stabilize by rotating them clockwise. Repeat the above steps for all wooden legs for height adjustment.

Backrest Adjustment

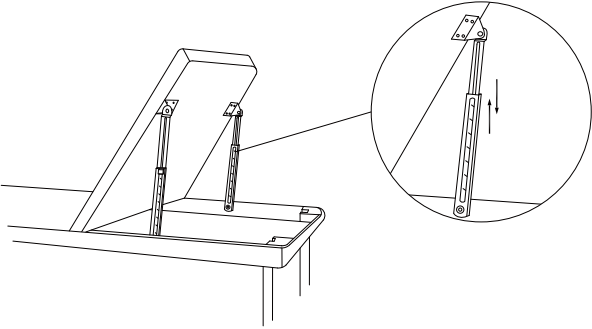

Lifting: Gently lift the backrest with 2 hands, and the lifter would respond accordingly, which allows the selection of desired height.

Putting back: Lift the backrest to the maximum height, and then put it down naturally.

- EN -

Product Maintenance

- a. The product is designed for massage, not for intense activities.
- b. Please wipe the surface of the product using mildly damp cloth or sponge for cleaning.
- c. Please do not clean the surface of the product with any corrosive detergent.
- d. Please fold and store the product if it's not to be used for long.

Warranty & Customer Support

If you need help regarding product issues, after-sales service, or other relevant problems, please contact us via the email on the cover.

* Note: Products purchased directly from Naipo benefit from our after-sales service. If the purchases are made from other selling teams, we recommend you contact them for warranty issues.

- EN -